

PRECISION SPINE
REFORM[®]
MODULAR PEDICLE SCREW

VERSATILITY & VISIBILITY

The Reform[®] Modular Pedicle Screw System provides surgeons increased flexibility, versatility and visibility to meet the varying requirements of degenerative and trauma procedures. Reform features a cobalt chrome tulip, a titanium triple lead proximally tapered thread to deliver strength, stability and efficiency to all thoracolumbar constructs. Modular tulips are available in standard and reduction configurations. Testing shows that the Reform Modular intraoperatively assembled construct has a pull off strength as strong as a Reform preassembled screw providing optimal security.




COBALT CHROME TULIP



REDUCTION TULIPS



TRIPLE LEAD THREAD

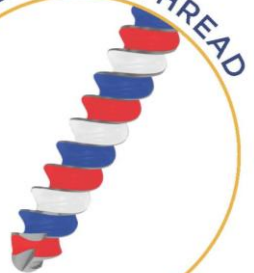

REFORM[®] MODULAR PEDICLE SCREW

Modular Screw Driver offers a secure, low profile interface optimizing maneuverability, control, and dependability.



Semi Aggressive Bone Planar approximates the outer diameter of the tulip to allow for bone decortication to promote fusion and proper seating of the tulip

Low Profile Tulip Clamp provides a stable and secure interface and is available in straight and angled configurations. An audible click ensures that the tulip is securely assembled

Polyaxial Screws provide 30° angulation in all planes for intraoperative flexibility.

Proximal Tapered Thread design increases Bone Screw interface, enhancing pull-out strength.


Triple Lead Thread allows for efficient screw delivery.

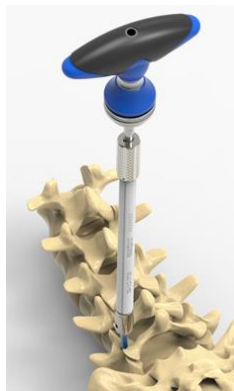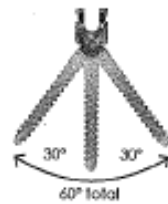

Screw Shank Driver



Bone Planar

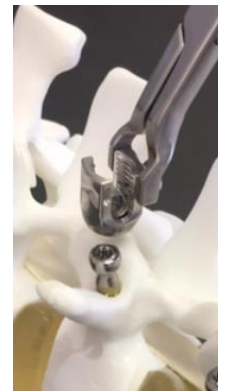

Tulip Clamp

Cobalt Chrome Tulip reduces risk of head splay and increases overall construct stabilization. Standard and reduction tulip options are available.



Square Threaded Locking Cap geometry reduces risk of cross threading

Hexalobular Drive Fittings on Bone Screw and Locking Cap reduces incidence of toggle and stripping



Reform Modular Sizing

4.5mm Diameter	5.5mm Diameter	6.5mm Diameter	7.5mm Diameter	8.5mm Diameter	9.5mm Diameter
25mm	30mm	30mm	35mm	40mm	60mm
30mm	35mm	35mm	40mm	45mm	70mm
35mm	40mm	40mm	45mm	50mm	80mm
40mm	45mm	45mm	50mm	55mm	
45mm	50mm	50mm	55mm	60mm	
		55mm		70mm	
				80mm	

Reform Rod Sizing – 5.5mm ø

Ti Straight	Ti Curved	CoCr Straight	CoCr Curved
80mm	35-80mm (5mm)	80mm	35-80mm (5mm)
100mm	90-120mm (10mm)	100mm	90-120mm (10mm)
120mm		120mm	
200mm		200mm	
300mm		300mm	
400mm		400mm	
500mm		500mm	

Precision Spine, Inc.
 2050 Executive Drive, Pearl, MS 39208
 Customer Service: 1.888.241.4773
 Phone: 601.420.4244
 Toll Free: 877.780.4370
 Fax: 601.420.5501
www.precisionspineinc.com



CE 2797