

PRECISION SPINE  
**REFORM**<sup>®</sup>  
PEDICLE SCREW SYSTEM

**STRENGTH, STABILITY AND EFFICIENCY**

Reform<sup>®</sup> is a comprehensive pedicle screw system that is designed to meet the varying requirements of degenerative, trauma and deformity procedures. Reform features a cobalt chrome tulip, a titanium triple lead proximally tapered thread and titanium and cobalt chrome rods to deliver strength, stability and efficiency to all thoracolumbar constructs. Reduction and uniplanar screw options, along with a full line of hooks, dominoes and offsets, complete the system to simplify the procedure and accommodate individual patient anatomy.



PRECISION SPINE<sup>®</sup>

*Discover the Difference*



# REFORM® PEDICLE SCREW SYSTEM



**Cobalt Chrome Tulip**  
reduces risk of head splay  
with DVR maneuvers

**Square Threaded  
Locking Cap**  
geometry reduces risk  
of cross threading



**Hexalobular Drive Fittings**  
on Bone Screw and Locking  
Cap reduces incidence of  
toggle and stripping




**Polyaxial Screws**  
provide 30° angulation in all  
planes for intraoperative flexibility

**Proximal Tapered Thread**  
design increases bone/screw  
interface enhancing pull-out  
strength



**Reduction Tower  
Quick Connect**  
allows for slip-on  
insertion over  
the tulip

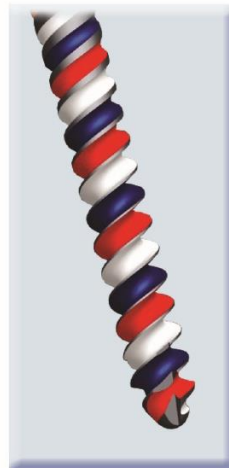

**Triple Lead Thread**  
allows for efficient  
screw delivery

## REFORM® SIZING

### Screw Sizes

4.5	25-45mm (5mm)
5.5	30-50mm (5mm)
6.5	30-55mm (5mm)
7.5	35-55mm (5mm)
8.5	40-60mm (5mm) 70 & 80mm
9.5	60-80mm (10mm)

### Rod Diameter

5.5mm

### Ti Straight

80, 100, 120, 200, 300, 400  
& 500mm

### Ti Curved

35-80mm (5mm)  
90-120mm (10mm)

### CoCr Straight

80, 100, 120, 200, 300, 400  
& 500mm

### CoCr Curved

35-80mm (5mm)  
90-120mm (10mm)

## Precision Spine, Inc.

2050 Executive Drive, Pearl, MS 39208  
Customer Service: 1.888.241.4773  
Phone: 601.420.4244  
Toll Free: 877.780.4370  
Fax: 601.420.5501  
www.precisionspineinc.com



CE 2797